



LILY Spinal Fluid Manometer

✓ Disposable ✓ Reliable ✓ Stable

Lily designs the spinal manometer as a tool connecting with lumbar puncture needle to provide a safe, reliable, flexible CSF pressure measurement. Product can be applied in emergency department, medical wards, and neurology centers.



LILY MEDICAL

A CSF analysis may include tests to diagnose :


- Infectious diseases of the brain and spinal cord : Meningitis and encephalitis
- Autoimmune disorders, such as Guillain - Barre Syndrome and Multiple Sclerosis (MS)
- Bleeding in the brain
- Brain Tumors

Swivel male lure adaptor on the 3-way stopcock provides a reliable connectivity during measuring.



CVP-S002
 Scale range : 4-45cmH2O
 Manometer length : 450mm
 Min. graduation : 2mm
 product content volume: <5ml
 4-Way stopcock

Transparent measuring rod allows easy and immediate inspection of CSF pressure during operation.

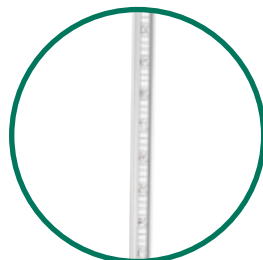

CVP-S001
 Scale range : 4-45cmH2O
 Manometer length : 450mm
 Min. graduation : 2mm
 product content volume : <5ml
 3-Way stopcock

Features

- The manometer tube has graduations marked from 4 cm to 45 cm and is attached to a stopcock.
- Disposable and for single patient use, reducing the risk of cross-infection in hospitals.
- Made of polycarbonate (PC) instead of glass for safer operation.
- Sturdy and one-piece construction, ready for immediate use.
- DEHP-free and latex-free; safe for patients.
- Passed hemolysis testing.

Order Code

Model Number	Description	Packing
CVP-S001	Spinal Manometer with 4-Ways stopcock.	Sterile, for single use, packed in cartons 25 units.
CVP-S002	Spinal Manometer with 3-Ways stopcock.	



SGS Certificate



LILY MEDICAL

TEL +886-37-612-625
 FAX +886-37-612-720
 E-MAIL customer.service@lily-medical.com
 WEB www.lily-medical.com
 ADD NO.28-2, Shun Jeau Diann, Chunan Town.
 Miaoli County, 35056. Taiwan R.O.C

